

Change Notice No.	Notification date	Implementation date
CN/F/21/007	December 31, 2021	January 31, 2022

Subject : Improved design of Cartridge Filters with adapter code 'Seal-M'

Scope : This change notification will affect the products with following catalog numbers: (# can be any alphabet or numeral):

- CB#####JO#####
- CP#####JO#####
- CW#####JO#####

Background:

In our efforts to provide the user, products with better functionality, certain changes have been made in the design of **mdj** Cartridge Filters with adapter code 'Seal-M' and having catalog numbers as mentioned above.

Need for Change:


The design of **mdj** Cartridge Filters with adapter code 'Seal-M' include a collar, at bottom of the cartridge element and just above the adapter used for the fitment of cartridge filter in respective housing. It was observed that in some of the cartridge housings, the collar of the fitted cartridge filter (with adapter code 'Seal-M') hinders the flow path by covering some part of fluid inlet of the housing. To overcome this issue, the design of **mdj** Cartridge Filters with adapter code 'Seal-M' has been changed and the outer diameter of collar is reduced.

No change has been done in any other critical dimensions or materials of construction or manufacturing process of **mdj** Cartridge Filters with adapter code 'Seal-M'.

Change Notice No.	Notification date	Implementation date
CN/F/21/007	December 31, 2021	January 31, 2022

How Does It Affect the User:

From a practical point of view the changes done in mdi Cartridge Filters with adapter code 'Seal-M' will enhance the functionality during actual use. All performance specifications and product attributes including the fitment with respective housing will remain unchanged, except the reduced diameter of collar of mdi Cartridge Filters with adapter code 'Seal-M'.

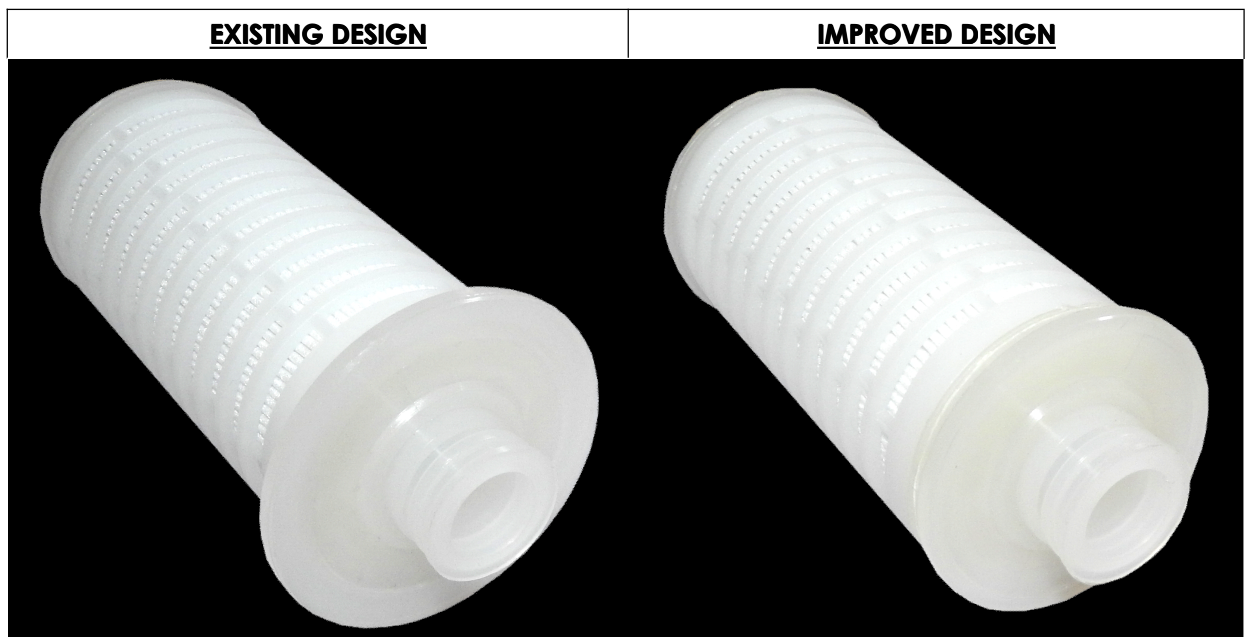

Implementation of Change:

mdi Cartridge Filters with adapter code 'Seal-M' and having catalog numbers as mentioned above, manufactured after the implementation date will have the said changes. However, the available stock of mdi Cartridge Filters with adapter code 'Seal-M' having old design and catalog numbers as mentioned above will be received by you till the stocks last.

In case you have any queries, please feel free to contact our Technical Support Team at 'info@mdimembrane.com'.

**Head - Quality Assurance
Advanced Microdevices Pvt. Ltd.**