

Change Notice No.	Notification date	Implementation date
CN/19/001	February 05, 2019	March 05, 2019

Subject : Change in the inlet/outlet connection of sampling bags used in mdi Single Use Systems

Scope : This change notification will affect the products with the following catalog number:

- A001XXXXXXXXSC301
- A002XXXXXXXXSC301
- A012XXXXXXXXXC301
- A028XXXXXXXXSS301
- A029XXXXXXXXSS301
- A039XXXXXXXXXC301
- A044XXXXXXXXSS301
- A075XXXXXXXXXC301
- A076XXXXXXXXXC301
- A077XXXXXXXXXC301
- A078XXXXXXXXXC301
- A138XXXXXXXXSC301
- A202XXXXXXXXXC301
- A203XXXXXXXXXC301
- A216XXXXXXXXXC301
- ABDXAMLXXXXX302
- ABDXBMLXXXXX302

Background:

In our efforts to provide the user more robust products, it has been decided to make certain changes in the design of the inlet/outlet connection of sampling bags used in **mdi** Single Use Systems with catalog numbers as mentioned above.


Need for Change:

Existing Sampling bags used in the **mdi** Single Use Systems with catalog numbers as mentioned above are manufactured using boat type inlet/outlet connection. These inlet/outlet connections are validated to withstand a pressure upto 5 psi. During usage of the Single Use Systems at customer end, sometimes the pressure generated on these sampling bags is more than the 5 psi. In such cases, a pathway is created next to the inlet/outlet connection of the sampling bag and leakage is observed through this pathway. To overcome this design constraint, it has been decided to introduce a new improved design of inlet/outlet connection in sampling bags. The new design of inlet/outlet connection has been verified for critical product performance parameters including increased operating pressure. No change has been done in any other dimensions or materials of construction.

Change Notice No.	Notification date	Implementation date
CN/19/001	February 05, 2019	March 05, 2019

How Does It Affect the User:

From a practical point of view this change only improves the operating pressure rating of the sampling bags used in Single Use Systems. All other performance specifications and product attributes remain unchanged, except the increased operating pressure.

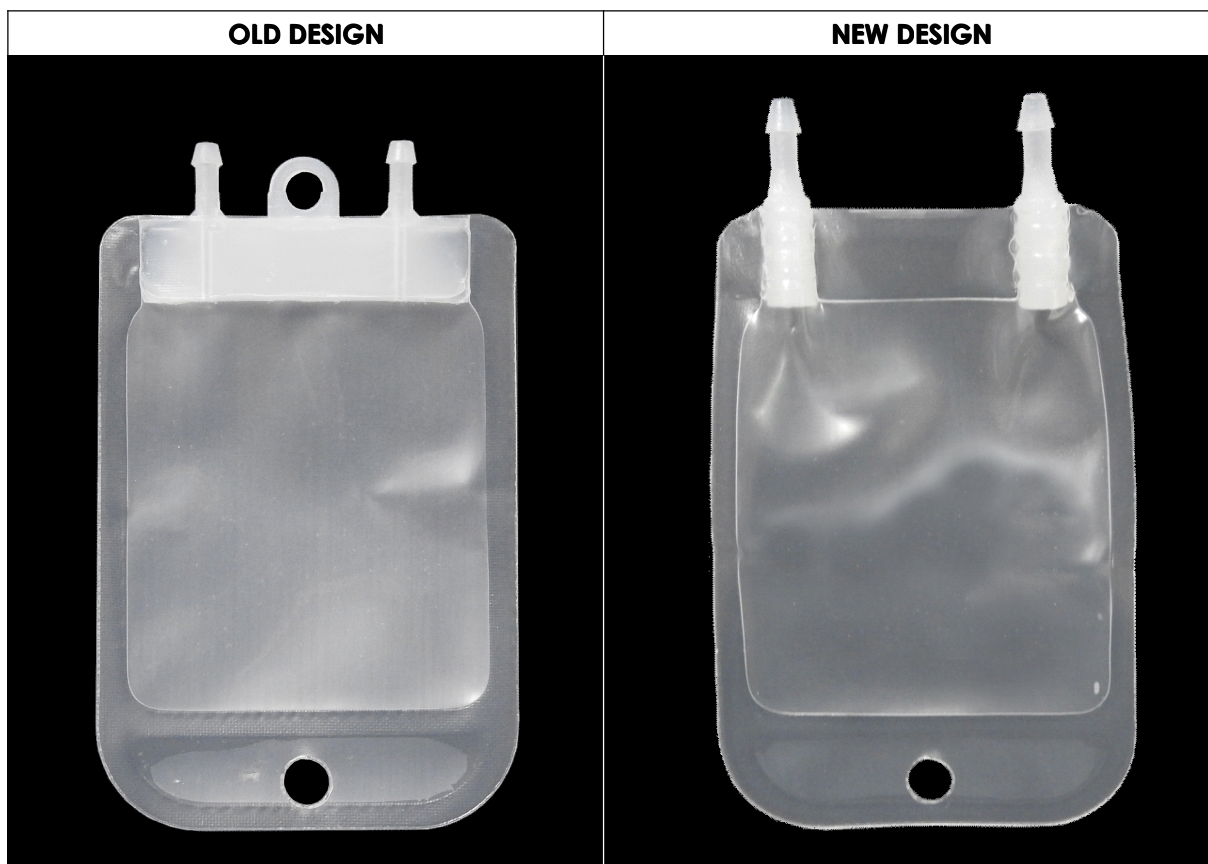

Effective Date:

mdi Single Use System with catalog numbers as mentioned above, manufactured after the implementation date (**March 05, 2019**) will have the said changes. However, the available stock of mdi Single Use System (with mentioned catalog numbers) having sampling bags with old type of inlet/outlet connection will be received by you till the stocks last.

In case you have any queries, please feel free to contact our Technical Support Team at info@mdimembrane.com



Quality Assurance
Advanced Microdevices Pvt. Ltd.